



**Corbac™
Adjustable
Contour
Back
Support**



**Forward™
Cushion**




**Control,
Comfort,
Protection**

It's all about contour...

Whether at work, rest, or play, the Ride Forward Cushion and Corbac surround the pelvis with supportive contours while reducing pressure around at-risk bony prominences. This helps improve postural control without compromise of skin protection.



Ride Designs® Forward™ Cushion

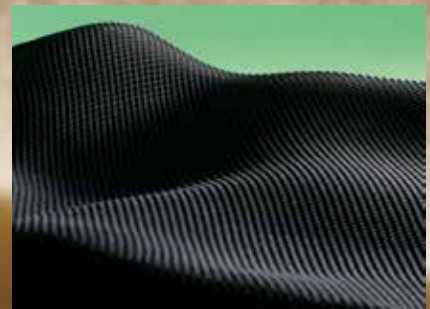

Control

Ride's patented contoured base helps promote and control improved posture and balance.

Comfort and protection

Top visco elastic layer provides protection and comfort as it conforms to body shape decreasing peak pressures.

Even more protection

Even under significant load, the spacer fabric outer cover does not compress. This allows air to move between the sitter and the cushion, helping to keep skin dry and cool.

Pediatric sizes available

Forward Cushions as small as 10 x 10" with a contour specially optimized for the pediatric pelvis.




▲ Spacer fabric in cover.



▲ Inner incontinent resistant cover, shown from the rear.

The cushion contour, spacer fabric outer cover and inner incontinent resistant cover all work together to keep the user and cushion dry. The inner cover really works. There are no seams or zippers in the path of liquid drainage, ensuring the cushion will not absorb moisture. Large accidents or spills simply pour off the back of the cushion.

Ride Designs® Corbac™ Adjustable Contour Back Support

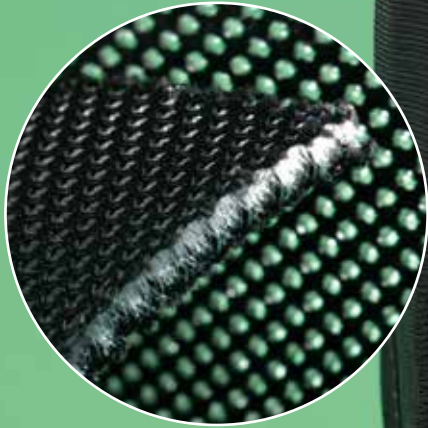

The unique spacer fabric material is used in both the Corbac and Forward Cushion for improved comfort and breathability.



▲ First, dial in the exact posterior contour with the Corbac's adjustment straps.

Then, refine the fit even further with the thoraco-lumbar-pelvic (TLP) supports. TLP supports provide precise adjustment of posterior and lateral contours. ▶

Lightweight

At two pounds, the Corbac™ weighs little more than standard back upholstery.

Precision fit

The Corbac is easily adjusted to help improve postural control and alignment. It does not decrease the wheelchair's depth or width. The sitter experiences fitting in, not sitting on, the wheelchair.



Reduce heat and moisture buildup

The spacer fabric overlay promotes air circulation and wicks moisture from the back surface.

Versatile

No screws, nuts, or tricky hardware, just a quick and easy mounting system that takes seconds to assemble onto the wheelchair.

A cinch to clean

The overlay detaches easily and can be machine washed.

Specs



Use CAM Wedges **symmetrically** to promote more upright posture and decrease ischial pressure.



Use CAM Wedges **asymmetrically** to help level a flexible pelvic obliquity.

Forward Cushion Specifications

Cushion Size (width x depth)		Part Number	Fits Wheelchair Widths		Fits Wheelchair Depth	
10 x 10"	25 x 25cm	FC-1010	9/10"	23/25cm	9/10"	23/25cm
10 x 12"	25 x 30cm	FC-1012	9/10"	23/25cm	11/12"	28/30cm
12 x 12"	30 x 30cm	FC-1212	11/12"	28/30cm	11/12"	28/30cm
12 x 14"	30 x 36cm	FC-1214	11/12"	28/30cm	13/14"	33/36cm
12 x 16"	30 x 41cm	FC-1216	11/12"	28/30cm	15/16"	38/41cm
14 x 14"	36 x 36cm	FC-1414	13/14"	33/36cm	13/14"	33/36cm
14 x 16"	36 x 41cm	FC-1416	13/14"	33/36cm	15/16"	38/41cm
16 x 16"	41 x 41cm	FC-1616	15/16"	38/41cm	15/16"	38/41cm
16 x 18"	41 x 46cm	FC-1618	15/16"	38/41cm	17/18"	43/46cm
16 x 20"	41 x 51cm	FC-1620	15/16"	38/41cm	19/20"	48/51cm
18 x 16"	46 x 41cm	FC-1816	17/18"	43/46cm	15/16"	38/41cm
18 x 18"	46 x 46cm	FC-1818	17/18"	43/46cm	17/18"	43/46cm
18 x 20"	46 x 51cm	FC-1820	17/18"	43/46cm	19/20"	48/51cm
20 x 16"	51 x 41cm	FC-2016	19/20"	48/51cm	15/16"	38/41cm
20 x 18"	51 x 46cm	FC-2018	19/20"	48/51cm	17/18"	43/46cm
20 x 20"	51 x 51cm	FC-2020	19/20"	48/51cm	19/20"	48/51cm

Ⓢ This size features pediatric cushion contour.

Ride CAM® Wedge Kit for Forward Cushion (FCD-WK)

Available for every cushion size.

Ride Corbac Sizes



Fits Wheelchair Widths		Corbac Height		Part Number
12"	30cm	12 - 14"	30 - 36cm	CB-1200
14"	36cm	12 - 15"	30 - 38cm	CB-1400
15"	38cm	13 - 16"	33 - 41cm	CB-1500
16"	41cm	14 - 17"	36 - 43cm	CB-1600
16"	41cm	17 - 20"	43 - 51cm	CB-1600T
17"	43cm	14 - 17"	36 - 43cm	CB-1700
18"	46cm	14 - 17"	36 - 43cm	CB-1800
18"	46cm	17 - 20"	43 - 51cm	CB-1800T
19 - 20"	48 - 51cm	14 - 17"	36 - 43cm	CB-1920
19 - 20"	48 - 51cm	17 - 20"	43 - 51cm	CB-1920T

▶ Note: For people at higher risk for skin breakdown, and/or more challenging postural control needs, consider the **Ride Custom Cushion and Back**, the ultimate in postural support and skin care.



Ride Designs®
a branch of Aspen Seating, LLC

phone: 303.781.1633
toll-free: 866.781.1633
www.ridedesigns.com

Patents: www.ridedesigns.com/patents
© 2015, Ride Designs. 090-033-D